



WRAPPEL® REUSABLE SURGICAL WRAPPERS

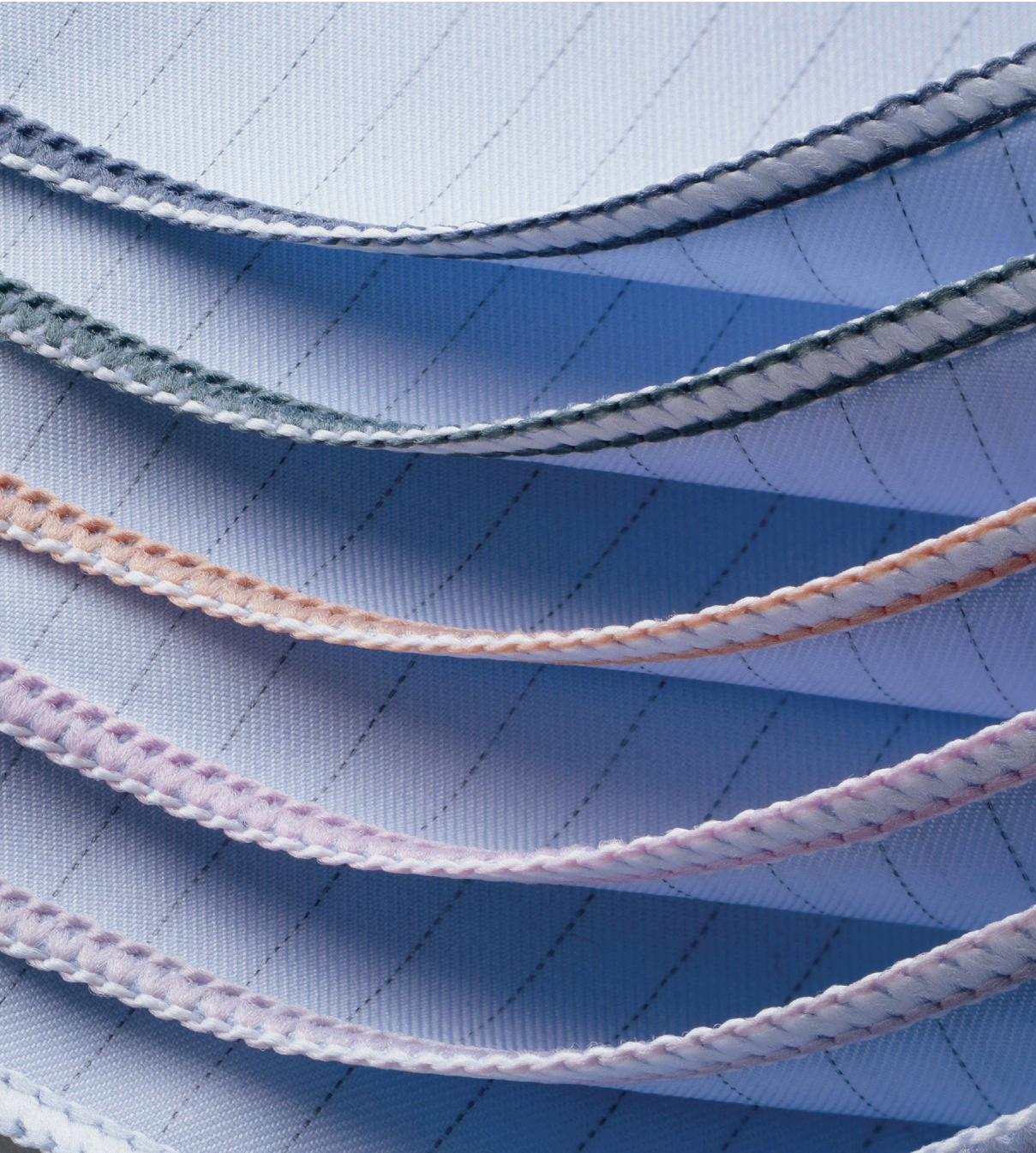

WRAPPEL®
Reusable Surgical Wrappers by Standard Textile

STANDARD
TEXTILE est. 1940

Reusable, dual-ply and liquid-resistant solution for wrapping surgical instrument sets, instrument trays and textile packs.

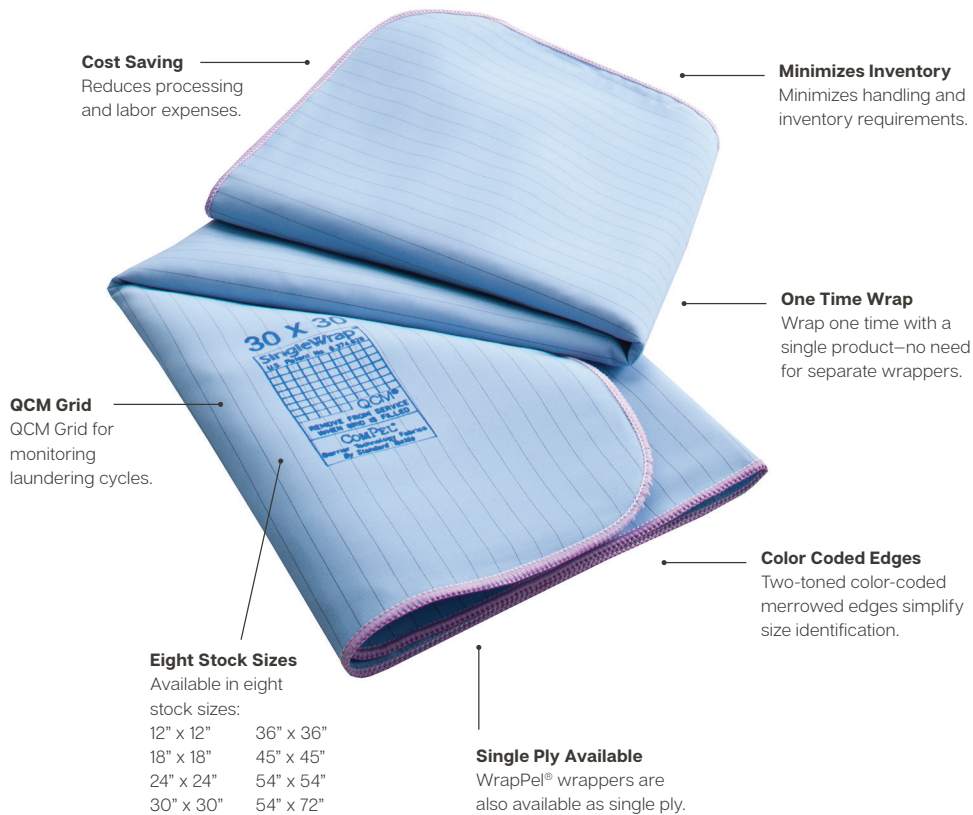

THE NEW STANDARD IN WRAPPING SURGICAL PACKS.

Patented WrapPel® SingleWrap® surgical wrappers are the industry's only dual-ply reusable, liquid-resistant surgical wrappers. Folding, inspection, and pack-making are simple and fast, resulting in substantial overall savings in labor, as well as reduced acquisition costs.

The multi-filament synthetic texturized yarns are virtually lint-free and the ripstop-weave fabric includes ESD (electro-static dissipative) yarns.

Standard Textile's WrapPel® Surgical Wrappers are liquid-resistant and virtually lint-free. They feature two-toned color-coded edges to help simplify size identification and a QCM grid for monitoring laundering cycles.

BEST PRACTICES:

To maximize performance, cost effectiveness, and longevity of reusable wrappers it is recommended to:

- Follow all appropriate AAMI and regulatory guidelines and instructions included in Standard Textile's IFU for WrapPel® wrappers.
- When wrapping instrument sets and other metal items, always utilize corner guards to protect wrappers from damage due to sharp edges and corners.
- When securing the wrapper with tape, create a tab at the end of the tape by folding the end over on itself. This will make it easier to remove tape from the wrappers before sending the used wrapper to the laundry for washing and processing.
- If/When a wrapper is damaged and needs to be patched or repaired, be sure to follow AAMI guidelines for the inspection and patching/repair of reusable surgical textiles.

For more information about surgical wrappers, or to contact a dedicated healthcare consultant, visit standardtextile.com or call 1.800.999.0400.

©2021. Standard Textile Co., Inc. All rights reserved.

SurgiTex
SURGICAL TEXTILE SOLUTIONS